

custom walker adjustments

With your arm extended down,
handle should be at wrist height.

With your hand placed on
the hand-grip, your elbow
should be flexed 20°-30°.



CAUTION:

Be aware of the following:

- Electrical cords in way
- Liquids on floor
- Throw rugs
- Toys or debris
- Furniture
- Doors and door steps
- Gaps in floor

WARNING:

1. Use wheels on front legs only or on all four legs with appropriate rear wheel brakes.
2. Do not push one side of the walker ahead and then other side without lifting the walker first.
3. Use Drive accessories only. Use of other manufacturers' accessories can result in improper fit and risk fall and/or injury.
4. Do not use on stairs or escalators.

lifetime limited warranty

Your Drive brand product is warranted to be free of defects in materials and workmanship for the lifetime of the original consumer purchaser.

This device was built to exacting standards and carefully inspected prior to shipment. This Lifetime Limited Warranty is an expression of our confidence in the materials and workmanship of our products and our assurance to the consumer of years of dependable service.

This warranty does not cover device failure due to owner misuse or negligence, or normal wear and tear. The warranty does not extend to non-durable components, such as rubber accessories, casters, and grips, which are subject to normal wear and need periodic replacement.

If you have a question about your Drive device or this warranty, please contact an authorized Drive dealer.

driveTM

© 2013 Medical Depot, Inc. All rights reserved.
Drive is a trademark of Medical Depot, Inc.
Port Washington N.Y. 11050 USA Made in China

universal folding
walker, two button
with 5" wheels

item # 10253-1




EU Authorized Representative
Drive Medical LTD
Aspley Industrial Estate
Elland, West Yorkshire
United Kingdom HX5 8JP

www.drivemedical.com

driveTM

universal walker instructions

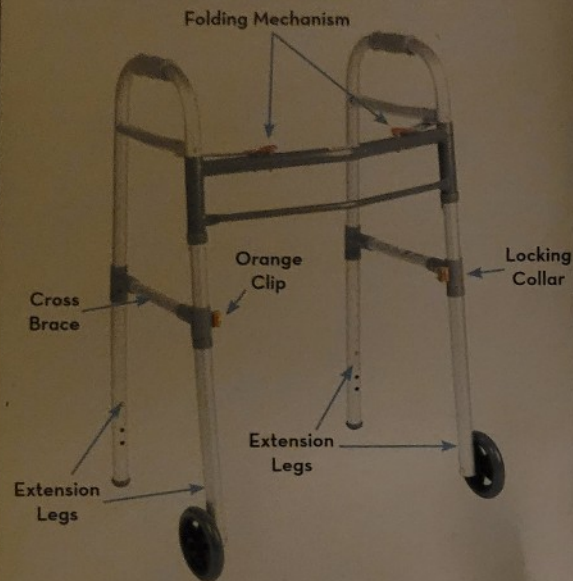

For use with all
Two Button, Universal Folding Walkers.

important set-up information

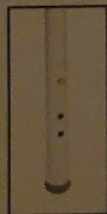


figure 1

1. Take out and unfold the Universal Deluxe 2-in-1 Folding Walker.
2. Adjust all four (4) extension legs (or wheels) to the required height (figure 1).



figure 2



figure 3

3. Detach locking collar by squeezing orange clip until the clip slides free of receptacle (figure 2 & 3).

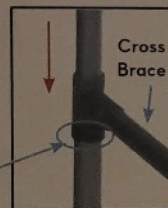


figure 4

4. Slide crossbrace to top of back leg extension collar (figure 4).

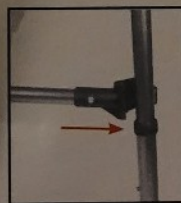


figure 5



figure 6

5. Pull black locking collar forward until flush with frame (figure 5). Close locking collar and insert orange clip into receptacle (figure 6).

NOTE: When adjusting the walker from a taller position to a lower position, first release the black locking collar with the orange clip. Slide the cross brace up, re-secure black locking collar with orange clip and set legs to desired height. The process may need to be repeated depending upon how low the adjustment is.

IMPORTANT: Test secure fit of orange clip by pulling on clip from front of walker.

wheelchair to walker transfer

Making certain wheels on the wheelchair are locked lift your body by pushing down on the armrests and carefully transfer one hand at a time to the walker.

wheelchair to walker transfer

(with lower extremity injury to right leg or foot)

Making certain wheels on the wheelchair are locked, lift your body by pushing down on the armrests and carefully transferring first the left hand, then the right hand to the walker. (If injury to left leg or foot transfer the right hand first, followed by left.)

three-point gait

(for walking with walker)

Step 1

Align the middle of your foot with the back legs of the walker.



Step 2

Slide the walker forward and place it comfortably in front of you.



Step 3

Carefully, walk forward one step at a time.

